

# **Export from GeneMapper ID-X & subsequent import into DNAXs/EuroForMix or Casesolver**

## **1. GeneMapper ID-X:**

- Check Samples and Genotypes Tabs
- If already available, choose Table Setting for DNAXs export or choose New... and create a corresponding template
- To check the setting of a table setting template go to Main Menu => Tools => Table Setting Editor or Shortcut (ctrl+T)
- Adapt Table Settings for Genotypes Tab to input format of various software tools (minimum 4 columns): Sample name, marker, allele 1, height 1, allele 2, height 2, [...]  
More details: Manual DNAXs 2.0, page 55, Sample Name Pattern, general file format description
- Main Menu => File => Export Table (= current view) OR Export Combined Table (combined table of samples and genotypes tab)

## **2. Import into DNAXs:**

Place GeneMapper export table (.txt or csv) in case folder within DNAXs

- In DNAXs import traces with either MUI Austria or INPS Lyon selected as 'sample pattern name' (for provided template including replicates) OR
- create a new 'sample pattern name' customized to your lab's sample ID method for import of your own Gene Mapper data

## **3. Import to EuroformMix:**

- Copy&paste relevant evidence profiles into a new text file to only include those into EuroForMix and any subsequent LR calculations (better overview, if not all evidence profiles are imported at once, but it is possible)

## **4. Import to Casesolver:**

- Delete allelic ladder, positive and negative profiles from Gene Mapper Textfile, then place text file in appropriate case folder
- Double-check: File => 'Set Working directory' (Fig. 1) AND Setup => 'set path with CASE folders' (Fig. 2) both include the same path (only then will you be able to see the import button)

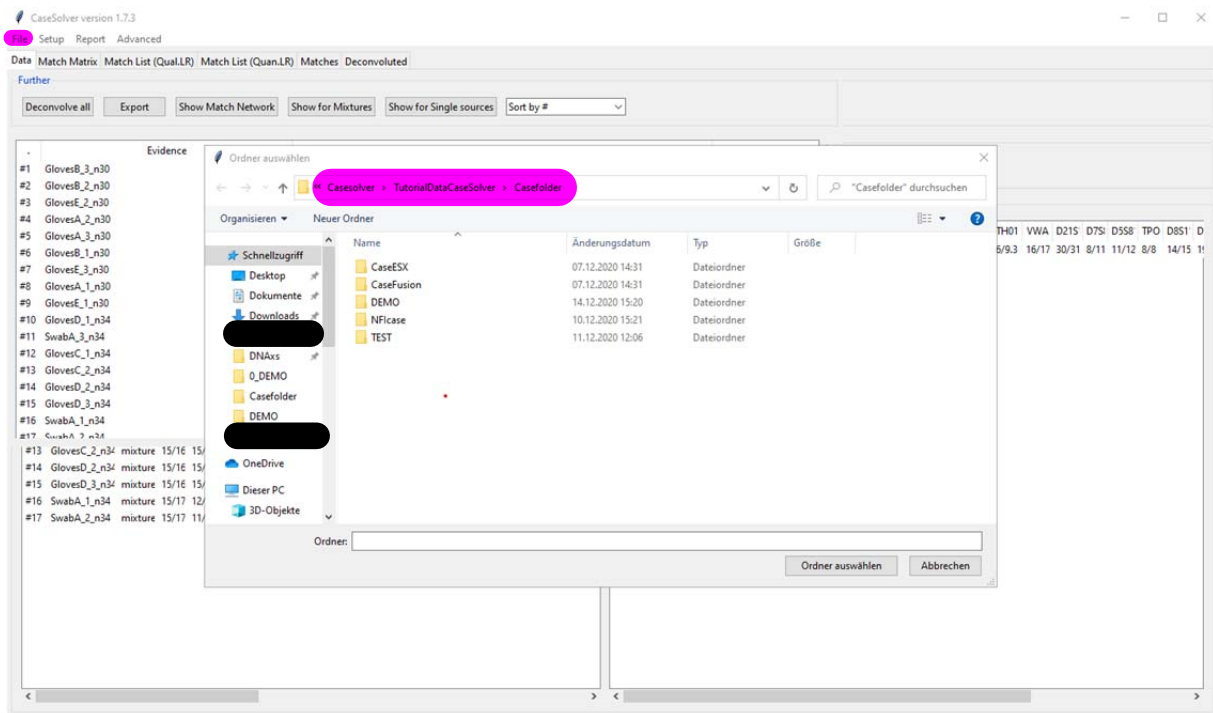


Figure 1: In Casesolver main menu 'File' => 'Set Working directory'

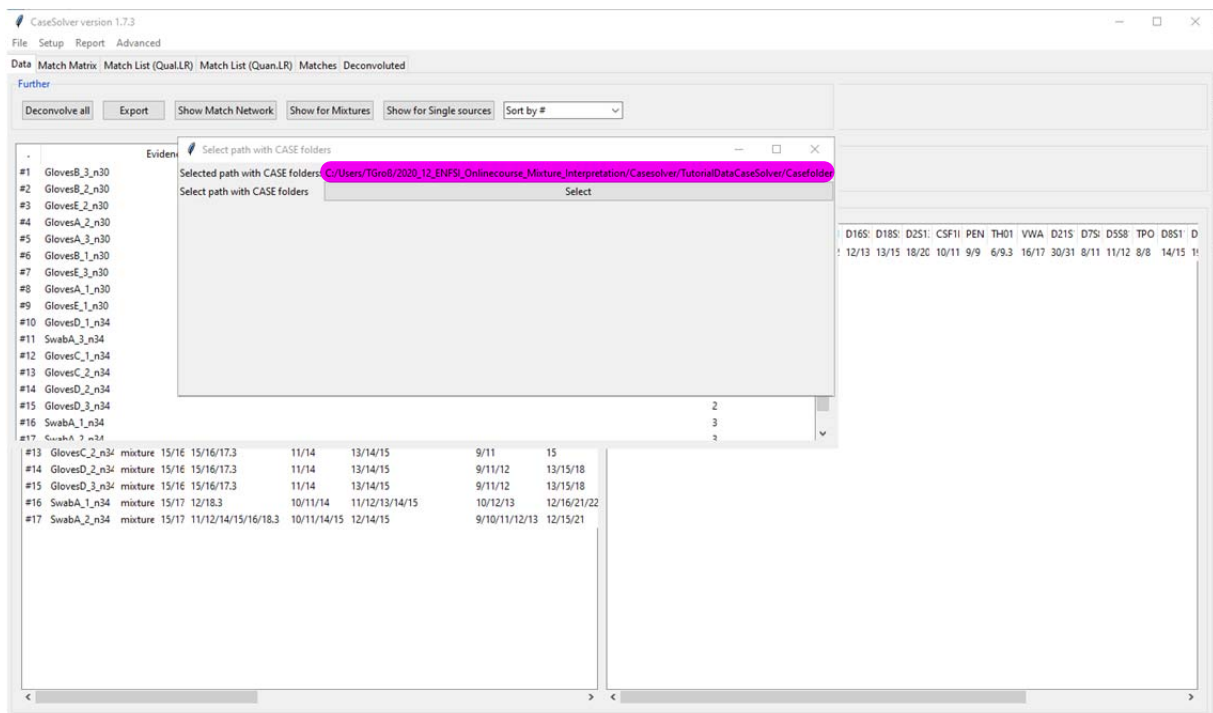


Figure 2: In Casesolver, Main Menu 'Setup' => 'set path with CASE folders'

	Status	Sample Name	Sample Type	Analysis Method	Panel	Size Standard
1		GlovesA_1_n30	Sample	3500_NGM_SE_2500_50bis100rfu	NGMSElect_panel_v1X	LIZ600_60-460
2		GlovesA_2_n30	Sample	3500_NGM_SE_2500_50bis100rfu	NGMSElect_panel_v1X	LIZ600_60-460
3		GlovesA_3_n30	Sample	3500_NGM_SE_2500_50bis100rfu	NGMSElect_panel_v1X	LIZ600_60-460
4		GlovesB_1_n30	Sample	3500_NGM_SE_2500_50bis100rfu	NGMSElect_panel_v1X	LIZ600_60-460
5		GlovesB_2_n30	Sample	3500_NGM_SE_2500_50bis100rfu	NGMSElect_panel_v1X	LIZ600_60-460
6		GlovesB_3_n30	Sample	3500_NGM_SE_2500_50bis100rfu	NGMSElect_panel_v1X	LIZ600_60-460
7		GlovesC_1_n34	Sample	3500_NGM_SE_2500_50bis100rfu	NGMSElect_panel_v1X	LIZ600_60-460
8		GlovesC_2_n34	Sample	3500_NGM_SE_2500_50bis100rfu	NGMSElect_panel_v1X	LIZ600_60-460
9		GlovesD_1_n34	Sample	3500_NGM_SE_2500_50bis100rfu	NGMSElect_panel_v1X	LIZ600_60-460
10		GlovesD_2_n34	Sample	3500_NGM_SE_2500_50bis100rfu	NGMSElect_panel_v1X	LIZ600_60-460
11		GlovesD_3_n34	Sample	3500_NGM_SE_2500_50bis100rfu	NGMSElect_panel_v1X	LIZ600_60-460
12		GlovesE_1_n30	Sample	3500_NGM_SE_2500_50bis100rfu	NGMSElect_panel_v1X	LIZ600_60-460
13		GlovesE_2_n30	Sample	3500_NGM_SE_2500_50bis100rfu	NGMSElect_panel_v1X	LIZ600_60-460
14		GlovesE_3_n30	Sample	3500_NGM_SE_2500_50bis100rfu	NGMSElect_panel_v1X	LIZ600_60-460
15		SwabA_1_n34	Sample	3500_NGM_SE_2500_50bis100rfu	NGMSElect_panel_v1X	LIZ600_60-460
16		SwabA_2_n34	Sample	3500_NGM_SE_2500_50bis100rfu	NGMSElect_panel_v1X	LIZ600_60-460
17		SwabA_3_n34	Sample	3500_NGM_SE_2500_50bis100rfu	NGMSElect_panel_v1X	LIZ600_60-460
18		Negativecontrol_1_n34	Negative Control	3500_NGM_SE_2500_50bis100rfu	NGMSElect_panel_v1X	LIZ600_60-460
19		NGMSELadder_1_n0	Allelic Ladder	3500_NGM_SE_2500_50bis100rfu	NGMSElect_panel_v1X	LIZ600_60-460
20		NGMSE_Ladder_1_n0	Allelic Ladder	3500_NGM_SE_2500_50bis100rfu	NGMSElect_panel_v1X	LIZ600_60-460
21		PositiveDNA007NGM-SE_1_n30	Positive Control	3500_NGM_SE_2500_50bis100rfu	NGMSElect_panel_v1X	LIZ600_60-460
22		PositiveDNA007NGM-SE_2_n30	Positive Control	3500_NGM_SE_2500_50bis100rfu	NGMSElect_panel_v1X	LIZ600_60-460

	Sample Name	Marker	Allele 1	Allele 2	Allele 3	Allele 4	Allele 5	Allele 6	Allele 7	Allele 8	Allele 9	Allele 10	Height 1	Height 2	Height 3	Height 4	Height 5	Height 6	Height 7	Height 8	Height 9	Height 10
1	GlovesA_1_n30	D10S1248	13	14	15								11580	1160	10586							
2	GlovesA_1_n30	vWA	15	17									6184	5139								
3	GlovesA_1_n30	D16S539	9	11	12								6552	206	3648							
4	GlovesA_1_n30	D2S1338	17	18	20	24							263	3795	3043	102						
5	GlovesA_1_n30	AMEL	X	Y									31944	238								
6	GlovesA_1_n30	D8S1179	13	14									21042	608								
7	GlovesA_1_n30	D21S11	30	30.2	31.2								8191	8023	234							
8	GlovesA_1_n30	D18S51	13	15									3985	3507								
9	GlovesA_1_n30	D22S1045	16	17									23277	1111								
10	GlovesA_1_n30	D19S433	13	14	14.2								871	13411	13414							
11	GlovesA_1_n30	TH01	6	7	8								9284	254	8953							
12	GlovesA_1_n30	FGA	19	21	22	23							146	598	8448	7182						
13	GlovesA_1_n30	D2S441	11	14									10055	11995								
14	GlovesA_1_n30	D3S1358	15	16									8625	8220								
15	GlovesA_1_n30	D1S1656	15	16	17.3								6975	7707	191							
16	GlovesA_1_n30	D12S391	20	22	23								4100	557	3807							
17	GlovesA_1_n30	SE33	19	22.2	28.2								2393	224	3242							
18	GlovesA_2_n30	D10S1248	13	14	15								13591	1353	12621							
19	GlovesA_2_n30	vWA	15	17									5702	5135								
20	GlovesA_2_n30	D16S539	9	11	12								5599	283	4781							
21	GlovesA_2_n30	D2S1338	17	18	20								571	4658	4436							
22	GlovesA_2_n30	AMEL	X	Y									32775	505								
23	GlovesA_2_n30	D8S1179	13	14									22104	451								
24	GlovesA_2_n30	D21S11	30	30.2	32.2								8100	7589	173							
25	GlovesA_2_n30	D18S51	13	15	18								4809	4654	193							
26	GlovesA_2_n30	D22S1045	16	17									29132	1339								
27	GlovesA_2_n30	D19S433	13	14	14.2								1264	13999	14748							
28	GlovesA_2_n30	TH01	6	7	8								11683	569	11068							
29	GlovesA_2_n30	FGA	21	22	23								536	8348	8321							
30	GlovesA_2_n30	D2S441	11	14									11573	11329								
31	GlovesA_2_n30	D3S1358	15	16									9971	8873								
32	GlovesA_2_n30	D1S1656	15	16	18								7700	7206	172							
33	GlovesA_2_n30	D12S391	16	20	22	23							119	4212	718	5069						
34	GlovesA_2_n30	SE33	19	28.2	30.2								2565	2822	158							
35	GlovesA_3_n30	D10S1248	13	14	15								17789	1527	15038							
36	GlovesA_3_n30	vWA	15	17									7110	8042								
37	GlovesA_3_n30	D16S539	9	11	12								5274	572	7064							
38	GlovesA_3_n30	D2S1338	17	18	20								301	3699	3238							
39	GlovesA_3_n30	AMEL	X	Y									32787	245								

	Sample Name	Marker	Allele 1	Allele 2	Allele 3	Allele 4	Allele 5	Allele 6	Allele 7	Allele 8	Allele 9	Allele 10	Height 1	Height 2	Height 3	Height 4	Height 5	Height 6	Height 7	Height 8	Height 9	Height 10
1	GlovesA_1_n30	D10S1248	13	14	15								11580	1160	10586							
2	GlovesA_1_n30	vWA	15	17									6184	5139								
3	GlovesA_1_n30	D16S539	9	11	12								6552	206	3648							
4	GlovesA_1_n30	D2S1338	17	18	20								263	3795	3043	102						
5	GlovesA_1_n30	AMEL	X	Y									31944	238								
6	GlovesA_1_n30	D8S1179	13	14									21042	608								
7	GlovesA_1_n30	D21S11	30	30.2	31.2								8191	8023	234							
8	GlovesA_1_n30	D18S51	13	15									3985	3507								
9	GlovesA_1_n30	D22S1045	16	17									23277	1111								
10	GlovesA_1_n30	D19S433	13	14	14.2								871	13411	13414							
11	GlovesA_1_n30	TH01	6	7	8								9284	254	8953							
12	GlovesA_1_n30	FGA	19	21	22	23							146	598	8448	7182						
13	GlovesA_1_n30	D2S441	11	14									10055	11995								
14	GlovesA_1_n30	D3S1358	15	16									8625	8220								
15	GlovesA_1_n30	D1S1656	15	16	17.3								6975	7707	191							
16	GlovesA_1_n30	D12S391	20	22	23								4100	557	3807							
17	GlovesA_1_n30	SE33	19	22.2	28.2								2393	224	3242							
18	GlovesA_2_n30	D10S1248	13	14	15								13591	1353	12621							
19	GlovesA_2_n30	vWA	15	17									5702	5135								
20	GlovesA_2_n30	D16S539	9	11	12								5599	283	4781							
21	GlovesA_2_n30	D2S1338	17	18	20								571	4658	4436							
22	GlovesA_2_n30	AMEL	X	Y									32775	505								
23	GlovesA_2_n30	D8S1179	13	14									22104	451								
24	GlovesA_2_n30	D21S11	30	30.2	32.2								8100	7589	173							
25	GlovesA_2_n30	D18S51	13	15	18								4809	4654	193							
26	GlovesA_2_n30	D22S1045	16	17									29132	1339								
27	GlovesA_2_n30	D19S433	13	14	14.2								1264	13999	14748							
28	GlovesA_2_n30	TH01	6	7	8								11683	569	11068							
29	GlovesA_2_n30	FGA	21	22	23								536	8348	8321							
30	GlovesA_2_n30	D2S441	11	14									11573	11329								
31	GlovesA_2_n30	D3S1358	15	16									9971	8873								
32	GlovesA_2_n30	D1S1656	15	16	18								7700	7206	172							
33	GlovesA_2_n30	D12S391	16	20	22	23							119	4212	718	5069						
34	GlovesA_2_n30	SE33	19	28.2	30.2								2565	2822	158							
35	GlovesA_3_n30	D10S1248	13	14	15								17789	1527	15038							
36	GlovesA_3_n30	vWA	15	17									7110	8042								
37	GlovesA_3_n30	D16S539	9	11	12								5274	572	7064							
38	GlovesA_3_n30	D2S1338	17	18	20								301	3699	3238							
39	GlovesA_3_n30	AMEL	X	Y									32787	245								

Genotypes Table Settings:

Column Settings:

	Show	Column	Filtering	Content
1	<input type="checkbox"/>	Sample File	Show All Records	
2	<input checked="" type="checkbox"/>	Sample Name	Show All Records	
3	<input type="checkbox"/>	Sample ID	Show All Records	
4	<input type="checkbox"/>	Run Name	Show All Records	
5	<input type="checkbox"/>	Panel	Show All Records	
6	<input checked="" type="checkbox"/>	Marker	Show All Records	
7	<input type="checkbox"/>	Dye	Show All Records	N/A
8	<input checked="" type="checkbox"/>	Allele	Show All Records	
9	<input type="checkbox"/>	Size	Show All Records	
10	<input checked="" type="checkbox"/>	Height	Show All Records	
11	<input type="checkbox"/>	Peak Area	Show All Records	
12	<input type="checkbox"/>	Data Point	Show All Records	
13	<input type="checkbox"/>	Mutation	Show All Records	
14	<input type="checkbox"/>	AE Reason For Change	Show All Records	
15	<input type="checkbox"/>	Marker Edit Comment (MEC)	Show All Records	
16	<input type="checkbox"/>	Allele Display Overflow (ADO)	Show All Records	N/A
17	<input type="checkbox"/>	Marker Edit (ME)	Show All Records	N/A
18	<input type="checkbox"/>	Off-scale (OS)	Show All Records	N/A

Show Hide Show All Hide All

Font Settings:

Font: Arial  
Size: 13

Sort by

Sample File  
☒ Ascending  
☐ Descending

Then by

Dye  
☒ Ascending  
☐ Descending

Then by

Marker  
☒ Ascending  
☐ Descending

Factory Defaults

Allele Settings

Number of alleles  
 Show number of alleles 10

Allele position  
☐ Keep Allele, Size, Height, Area, Data Point, Mutation and Comment together

New Project...  
Ctrl+N

Open Project...  
Ctrl+O

Save Project  
Ctrl+S

Save Project As...

Export Table...  
Ctrl+E

Export Table for CODIS...

Export Combined Table...

Page Setup...

Print...  
Ctrl+P

Print Preview...

Project Options

Log Out

Exit  
Alt+F4

Genotypes																				
	Marker	Allele 1	Allele 2	Allele 3	Allele 4	Allele 5	Allele 6	Allele 7	Allele 8	Allele 9	Allele 10	Height	Height	Height	Height	Height	Height	Height	Height	Height
	D2S1338	17	18	20								621	6964	6872						
	AMEL	X	Y									32759	227							
	D8S1179	13	14									26130	418							
	D21S11	28	30	30.2								193	6744	12297						
	D18S51	13	15									9720	10009							
	D22S1045	16	17									25420	1269							
	D19S433	13	14	14.2								1004	14065	15237						
	TH01	6	7	8								11802	344	13907						
	FGA	21	22	23								730	9280	10251						
	D2S441	11	14									13386	11217							
	D3S1358	15	16									11666	11256							
	D1S1656	15	16	18								8352	7695	335						
100	GlovesB_3_n30	D12S391	20	22	23							6457	667	5443						
101	GlovesB_3_n30	SE33	19	28.2								6984	6019							
102	GlovesB_3_n30	SE33	19	28.2								6984	6019							
103	GlovesC_1_n34	D10S1248	13	14	15							5786	5731	14643						
104	GlovesC_1_n34	vWA	15	17								8653	12452							
105	GlovesC_1_n34	D16S539	9	11	12							4553	2676	5575						
106	GlovesC_1_n34	D2S1338	17	18	20	24						7148	2825	3944	2737					
107	GlovesC_1_n34	AMEL	X	Y								9719	2929							
108	GlovesC_1_n34	D8S1179	13	14								9587	13114							
109	GlovesC_1_n34	D21S11	28	30	32.2							3768	8699	6967						
110	GlovesC_1_n34	D18S51	13	15	18							4553	7494	4586						
111	GlovesC_1_n34	D22S1045	16	17								25358	9413							
112	GlovesC_1_n34	D19S433	13	14	14.2							5879	9660	10910						
113	GlovesC_1_n34	TH01	6	7								3448	3625							
114	GlovesC_1_n34	FGA	19	21	22	23						2616	4404	7988	5776					
115	GlovesC_1_n34	D2S441	11	14								13839	5654							
116	GlovesC_1_n34	D3S1358	15	16								9153	4933							
117	GlovesC_1_n34	D1S1656	15	16	17.3							2245	2703	2569						
118	GlovesC_1_n34	D12S391	16	20	22	23						4724	997	1151	1863					
119	GlovesC_1_n34	SE33	19	22.2	28.2	30.2						1068	8272	683	1906					
120	GlovesC_2_n34	D10S1248	13	14	15							2483	3351	8222						
121	GlovesC_2_n34	vWA	15	17								4826	8855							
122	GlovesC_2_n34	D16S539	9	11								2396	3430							
123	GlovesC_2_n34	D2S1338	17	18	20	24						1006	1209	2317	2338					
124	GlovesC_2_n34	AMEL	X	Y								9019	123							
125	GlovesC_2_n34	D8S1179	13	14								2502	4906							
126	GlovesC_2_n34	D21S11	30	30.2	32.2	33.2						2252	544	6885	223					
127	GlovesC_2_n34	D18S51	15									3731								



New Project... Ctrl+N

Open Project... Ctrl+O

Save Project Ctrl+S

Save Project As...

Export Table... Ctrl+E

Export Table for CODIS...

Export Combined Table...

Page Setup...

Print... Ctrl+P

Print Preview...

Project Options

Log Out

Exit Alt+F4

Table Setting: DNAXs																					
Genotypes																					
Marker	Allele 1	Allele 2	Allele 3	Allele 4	Allele 5	Allele 6	Allele 7	Allele 8	Allele 9	Allele 10	Height 1	Height 2	Height 3	Height 4	Height 5	Height 6	Height 7	Height 8	Height 9	Height 10	
D10S124	13	14	15								11580	1160	10586								
vWA	15	17									6184	5139									
D16S539	9	11	12								6552	206	3648								
D2S1338	17	18	20	24							263	3795	3043	102							
AMEL	X	Y									31944	238									
D8S1179	13	14									21042	608									
D21S11	30	30.2	31.2								8191	8023	234								
D18S51	13	15									3985	3507									
D22S104	16	17									23277	1111									
D19S433	13	14	14.2								871	13411	13414								
TH01	6	7	8								9284	254	8953								
FGA	19	21	22	23							146	598	8448	7182							
D2S441	11	14									10055	11995									
D3S1358	15	16									8625	8220									
D1S1656	15	16	17.3								6975	7707	191								
D12S391	20	22	23								4100	557	3807								
SE33	19	22.2	28.2								2393	224	3242								
D10S124	13	14	15								13591	1353	12621								
vWA	15	17									5702	5135									
D16S539	9	11	12								5599	283	4781								
D2S1338	17	18	20								571	4658	4436								
AMEL	X	Y									32775	505									
D8S1179	13	14									22104	451									
D21S11	30	30.2	32.2								8100	7589	173								
D18S51	13	15	18								4809	4654	193								
D22S104	16	17									29132	1339									
D19S433	13	14	14.2								1264	13999	14748								
TH01	6	7	8								11683	569	11068								
FGA	21	22	23								536	8348	8321								
D2S441	11	14									11573	11329									
D3S1358	15	16									9971	8873									
D1S1656	15	16	18								7700	7206	172								
D12S391	16	20	22	23							119	4212	718	5069							
SE33	19	28.2	30.2								2565	2822	158								
D10S124	13	14	15								17789	1527	15038								
vWA	15	17									7110	8042									
D16S539	9	11	12								5274	572	7064								
D2S1338	17	18	20								301	3699	3238								
AMEL	X	Y									32787	245									